

REMEDIUUM
Sky
SERIES

New Standard of HA Fillers

© Forever 18 INTERNATIONAL
www.remediumkorea.com

 ReMediumkorea  ReMedium  Remedium Korea

Fax: +82-2-2088-1029
Phone: +82-70-7525-1021
E-mail: contact@remediumkorea.com
Address: #1001 Daegil Medical Building, 384 Sagajeong-ro,
Jungnang-gu, Seoul, South Korea



REMEDIUUM
Sky
SERIES

1.1ml Hyaluronic acid with lidocaine

© Forever 18 INTERNATIONAL

SKY'S SECRET: YOUR JOURNEY TO TIMELESS BEAUTY



OUR MISSION: A BEAUTY JOURNEY PARTNER

At ReMedium, we stand by the principle that
"YOUR BEAUTY SOARS SKY-HIGH NOW"

We are dedicated to delivering safe and effective face filler solutions that enhance your natural features while promoting self-assurance.

We believe that everyone should feel comfortable and beautiful in their own skin, at any age. Our commitment is to provide you with top-quality face filler products that can help you achieve the aesthetic results you desire. Your beauty, your choice – we're here to support your journey to a more confident and radiant you.

SIGNATURE OF GLOBAL AESTHETICS

QUALITY GUARANTEED

ReMediumSKY products are manufactured in line with strict international standards- **ISO 13485:2016** and No. 22-B-0568 Rev.0 (for medical devices). All of our products are subject to strict quality controls, individually inspected, manually and visually assessed, one by one, before their release to the market.

CONTINUOUS INNOVATION

ReMediumSKY is committed to continuous innovation and dedicates a significant percentage of its yearly revenues to **R&D**.

STUDIES & CLINICAL TRIALS

To ensure the safety and efficacy of all ReMediumSKY products, clinical trials and tests are continuously conducted. Furthermore, our **Forever18 International Co., Ltd** is dedicated to regularly train and educate all ReMedium Family physicians around the Globe on how to apply the precise protocols in order to guarantee the best results for patients.



REFRESH
REVIVE

REMEDIUMSKY

14
•
15



WHERE BEAUTY MEETS GRACE

Be Confident, Be You

Your face is a canvas of expression, revealing your story, confidence, and unique beauty. It's time to unlock a new chapter with ReMediumSKY.

Get ready to show more of you!

Discover a personalized transformation with ReMediumSKY face fillers, designed to refresh, revive, and reveal your true essence. Our treatments add a natural, youthful glow, restoring your facial harmony and rejuvenating your inner radiance.

Elevate your beauty and embrace a renewed you. Experience the ReMediumSKY difference, where your true self shines through.



THE NEW STANDARD OF HA FILLERS




HIGH QUALITY RAW MATERIALS

ReMediumSKY uses only high-quality and safe materials listed in EDQM (European Directorate for the Quality of Medicines & Healthcare) for HA, lidocaine, and BDDE.

PURIFICATION PROCESS (Minimized BDDE residues)

ReMediumSKY minimizes impurities such as incompletely cross-linked BDDE (pendant and unreacted BDDE) and endotoxin residues through a 14-day purification process.

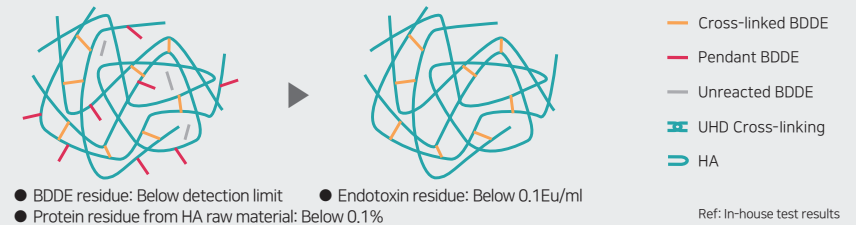

DISINTEGRATION TEST

After 48 hours of injecting hyaluronidase into 0.5ml of ReMediumSKY, it is completely disintegrated, proving its safety (Observed for 48 hours without shaking or mixing after the injection).

	0 hour	After 3 hours	After 6 hours	After 12 hours	After 48 hours
ReMediumSKY					

※ If side effects occur with hyaluronic acid fillers, use 200 IU or more of hyaluronidase and gently massage the treatment area. Please do not exercise too much pressure.

Ref: In-house research result

OPTIMIZED EASINESS IN USAGE

MANUFACTURING PROCESS

The differentiated manufacturing process increases the viscoelasticity of the filler and ensures safety.




25°C Assures safety by maintaining consistent temperature at 25°C for all manufacturing processes except sterilization;

14 days Minimizes residual BDDE and improves quality with 14 days purification process;

80% UHD Technology assures safety and high viscoelasticity with only 80% of BDDE compared to other companies.


UHD TECHNOLOGY (*KIPO Patent 10-2186646)

UHD (ULTRA HIGH DENSITY) is our patented Cross-linking Technology. It increases the viscoelasticity by improving the internal density at an optimal frequency for a certain period of time.

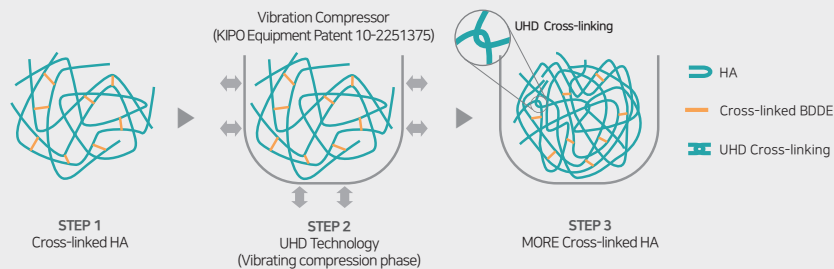

VISCOELASTICITY INCREASE RATE

The viscoelasticity (%) value of the product is increased using the UHD method.

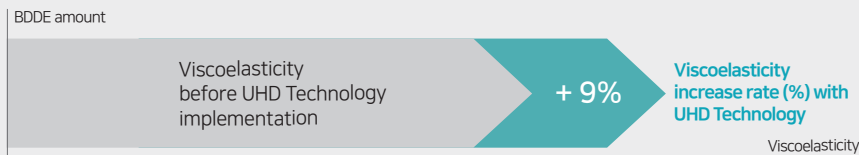

Ref: S.Thepharm Co.,Ltd, Method for producing a composition for filler using a vibrator for filler composition, KIPO Patent 10-2186646
Ref: KIPO Patent: (Patented by the Korean Intellectual Property Office)

Low extrusion force
Due to conical shape*

Low patient infection risk **
Higher product quality in terms of extractables and leachables contamination

Safe application
Extra grip surface at backstop and plunger rod for better hold



Comfortable handling
Solid designed components for a safe and pleasant routine work

Visual orientation
Easy and fast overview

Innovative design
Innovative user-centric approach lies at the core of our product design.


* EPC (End Product Control): Extrusion force is tested on every batch.

** Evaluation of component stability in accordance with international standards.

PROVEN EXCELLENCE IN EFFECTIVENESS

CLINICAL TRIAL RESULTS

ReMediumSKY clinical results show significant reduction of wrinkles compared to the control panel, with a high satisfaction level.

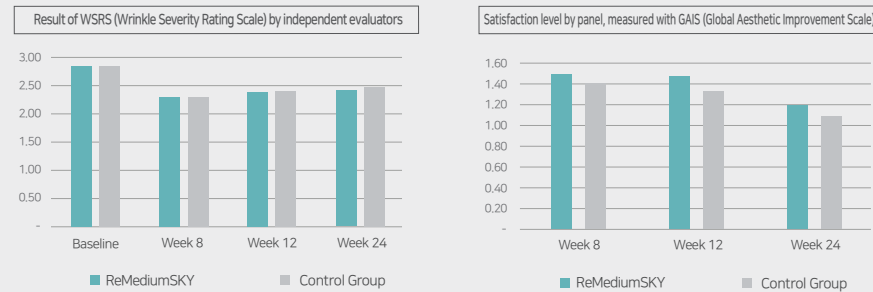


The result of the WSRS by an independent evaluator shows the degree of wrinkle depth improvement from the baseline at 8 weeks, 12 weeks, and 24 weeks.

The overall appearance improvement by the subject GAIS was found to be similar to or slightly higher than the control group.

Ref: Randomized multiple institutions to evaluate the efficacy and safety of ReMediumSKY Ultra and the control group for improvement of nasolabial folds and conducted clinical tests for double-blind (evaluator-subject), non-inferiority, counterpart, and confirmatory (Testing institution: Plastic Surgery at Severance Hospital, Dermatology at Severance Hospital, and Plastic Surgery at ASAN Medical Center).

EFFECTIVENESS

High cohesiveness of ReMediumSKY proves its effectiveness by pushing the FHL (Forehead horizontal lines) up to the level of surrounding tissues by taking the form of a small cylindrical filament. Residue remains after 12 months, ensuring long-term effectiveness.

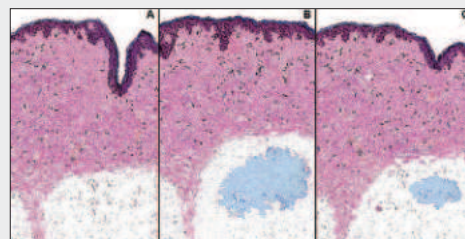


- A. Intradermal FHL depth before HA application stained with hematoxylin-eosin.
- B. No wrinkle in the dermis and small nodules formed in the subcutaneous fat, 48 hours after HA application.
- C. Partial recovery of the intradermal FHL depth and reduction in the nodule diameter in the subcutaneous fat, 12 months after application.

Ref: "Histological and Clinical Appraisal of Distribution and Long Acting New Hyaluronic Acid into the Subcutaneous Fat of the Forehead to Improve Horizontal Lines", ANTONIO CARLOS ABRAMO (BRAZIL).

DURABILITY

ReMediumSKY Ultra has proved its durability with an MRI scan. Despite the decrease in size and denseness of the spindle-shaped nodule, its presence in the subcutaneous fat after 12 months evidence the efficacy and long-acting of the HA filler.

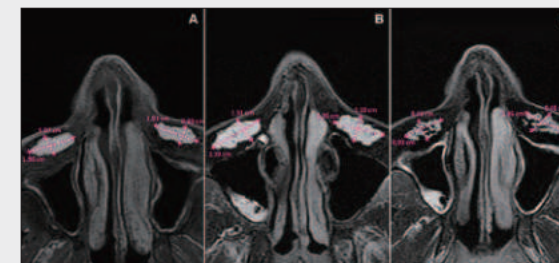


Measures of the largest axes of the spindle-shaped nodule

- A. At 24 hours.
- B. 1 month after HA application
- C. 12 months after HA application

Ref: "Effectiveness and Long-Acting of the Hyaluronic Acid Injected into the Nasolabial Fold Measured Through Magnetic Resonance Imaging", ANTONIO CARLOS ABRAMO (BRAZIL).

COHESIVITY

ReMediumSKY Ultra coheres better in subcutaneous fat than the dermis. When injected into the subcutaneous layer of the nasolabial fold, it functions as an effective and long-lasting volumizer attributable to its high degree of concentration and cohesivity.

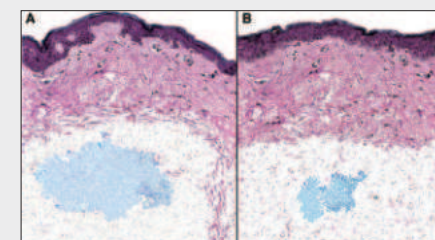

- A. A Nodule diameter 72 h after application into the subcutaneous fat.
- B. Moderate decrease in the nodule diameter 9 months after application

Ref: "Histological and Quantitative Appraisal of the Intradermal and Subcutaneous Distribution of the Hyaluronic Acid Injected into the Nasolabial Fold", ANTONIO CARLOS ABRAMO (BRAZIL).

FILLER INJECTION MAPPING



Effective Solutions

ReMediumSKY Deep

ReMediumSKY Soft

ReMediumSKY Soft

ReMediumSKY Deep

ReMediumSKY Deep


ReMediumSKY Soft

ReMediumSKY Ultra

ReMediumSKY Ultra

ReMediumSKY Ultra

ReMediumSKY Deep

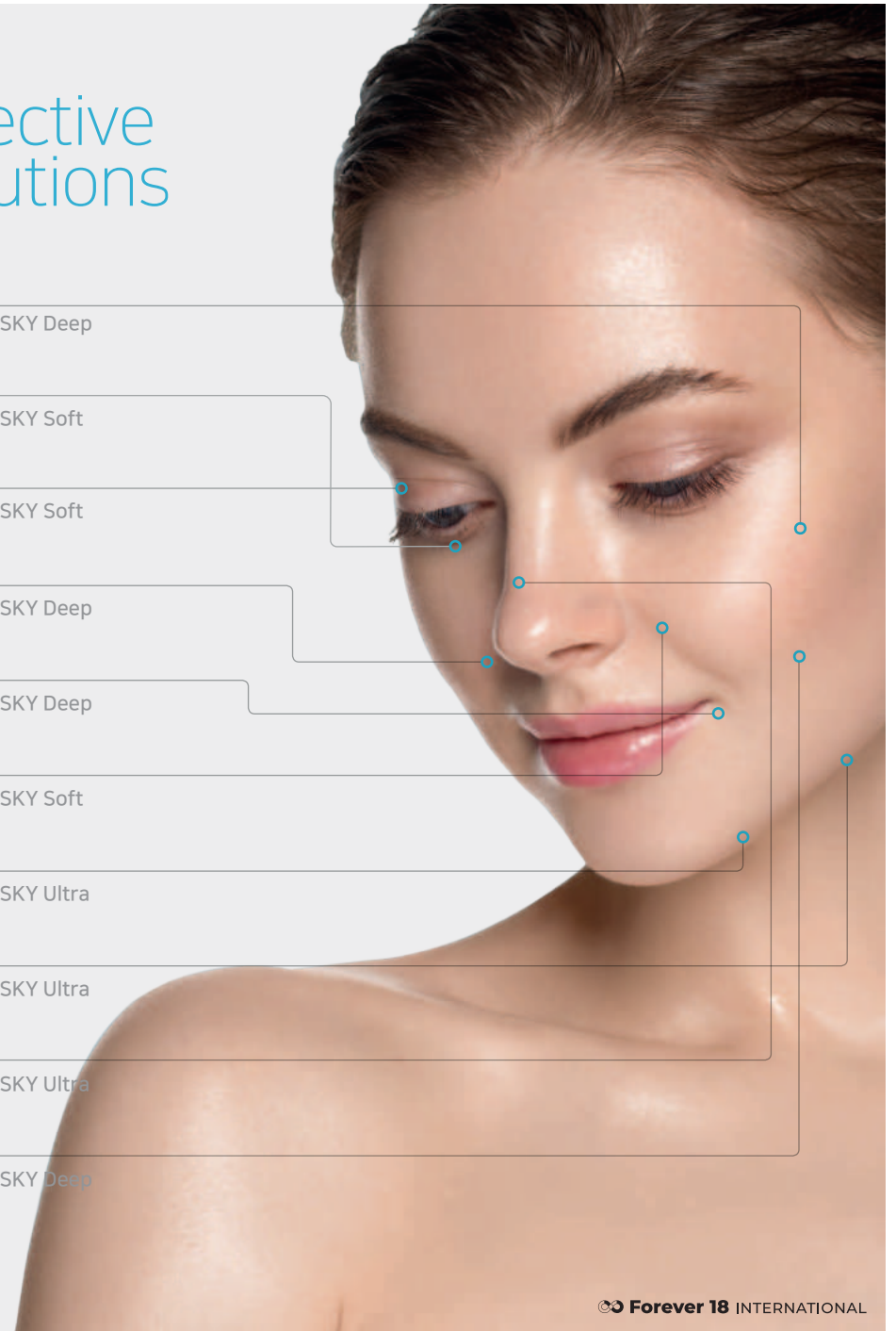




YOUR BEAUTY SOARS

SKY-HIGH NOW

Certified by EU, we produce Quality Monophasic Filler

TOTAL SOLUTION

Product Name	ReMediumSKY Soft	ReMediumSKY Deep	ReMediumSKY Ultra
HA concentration	24mg/ml	24mg/ml	24mg/ml
Needle Gauge	30G	27G	25G
Application & Indication	Lip defects or volume enhancement	Mid to deep dermis injection for correction of scars or deep wrinkles or facial tissue augmentation	A facial tissue augmentation for mid to deep dermal implantation for the correction of moderate to severe facial wrinkles and folds
Longevity	6 ~ 9 months	9 ~ 14 months	12 ~ 18 months
Elasticity modulus	50~120 pa	120~250 pa	Over 250 pa
Shape			

- The type of product is classified according to the difference in viscoelasticity.
- Osmotic pressure: 200~400mOsm/kg
- Endotoxins level: on average 0,1~0,5 EU/ml
- pH: 6,8~7,5

